




CAPITAL MEDICAL SUPPLIES PTY LTD
71 HEFFERNAN STREET
MITCHELL. ACT. 2911


PH. 02 62415511
FAX. 02 62417538
E-MAIL. info@capitalmedicalsupplies.com.au
ABN. 33 008 592 836

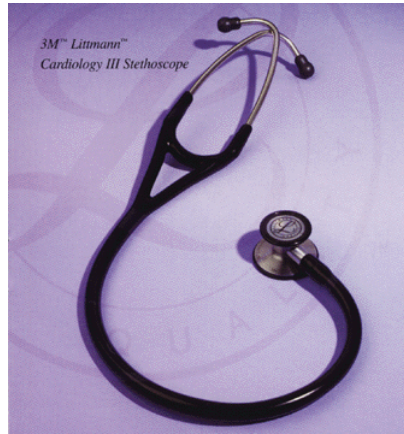


LITTMANN CARDIOLOGY 3 STETHOSCOPE

Product Information

The Littmann Cardiology III stethoscope is the very latest innovation in auscultation technology. Two patented tunable diaphragms allow for listening to both low and high frequency sounds with either side of the chestpiece. The large side can be used for auscultating adult patients while the small side is especially useful with pediatric or thin patients.

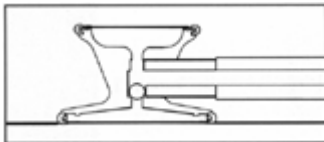

Bell Mode (low frequency)
Light contact with chestpiece allows you to hear low-frequency sounds.

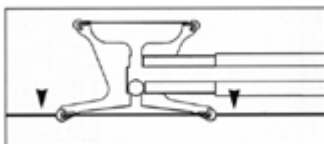

Diaphragm Mode (high frequency)
Without turning the chestpiece over, use firm pressure to hear high-frequency sounds.

Two Tubes-In-One Design

The Littmann Cardiology III tubing has been designed for acoustical excellence. Two tubes are incorporated into a single-tube design to eliminate noise artifact from tubes rubbing together.

Warranty

The quality and performance of the Littmann Cardiology III stethoscope is backed by a five-year warranty. After the warranty period, lifetime service and refurbishment is available.